

# SEPARATION REPORT

## Ultraspeed GPC Semi-micro Column TSKgel SuperHZ Series

### Table of Contents

1. Introduction	1
2. Features of TSKgel SuperHZ Series	1
3. Basic Properties	1
3-1. Resolution	2
3-2. Dependence of Height Equivalent to a Theoretical Plate (HETP) on Flow Rate	3
3-3. Effect of Sample Injection Volume	5
3-4. Effect of Sample Concentration	7
3-5. Optimization of Packing Material Particle Size	11
4. Application	12
5. Conclusion	14

## 1. Introduction

Ever since its emergence, gel permeation chromatography (GPC) has developed to a high-precision, high-reproducibility system with higher functions of the device, columns and packing materials. Tosoh Corporation introduced a group of ultra-high-speed organic solvent GPC columns, "TSKgel SuperH series," to the market ahead of the world, and realized a short-time, high-resolution, and solvent-saving GPC. Since we have commercialized a group of organic solvent semi-micro GPC columns, "TSKgel SuperHZ series" with higher performance in addition, this article introduces its basic features as well as application data.

## 2. Features of TSKgel SuperHZ Series

The SuperHZ series consists of SuperHZ1000 - 4000 in which a packing material of a different pore size is packed for each grade, and SuperH2M-N, M-M, and M-H in which packing materials with different pore sizes were mixed and the range of measurable molecular weights is extended (**Table-1**).

SuperHZ1000 - 4000 are capable of measuring from monomers, polymer additives, and oligomers to polymers up to molecular weight of several hundred thousands by selecting the grade. Since the resolution of the column is essential in this range of molecular weight (especially in the low molecular weight range), fine particles (3 $\mu$ m) are used for packing.

On the other hand, the mixed grades (SuperH2M-N, M-M, and M-H) are capable of measuring from oligomers to polymers with molecular weights up to tens of millions by selecting the grade. Packing materials of appropriate particle sizes are used for packing, in order to place focus on resolution in the low-molecular weight range and to avoid shear degradation of polymer in the high-molecular weight region. The features of SuperHZ series are shown in **Table-2**. In addition, calibration curves of SuperHZ series are shown in **Figures-1** and **-2**.

**Table-1 List of TSKgel SuperHZ Series**

Grade	Exclusion limit (polystyrene)	Particle size ( $\mu$ m)	Theoretical plates (TP/15cm)	Column size (mm I.D $\times$ cm)
TSKgel SuperHZ1000	$1 \times 10^3$	3	16,000	4.6 $\times$ 15
				6.0 $\times$ 15
TSKgel SuperHZ2000	$1 \times 10^4$	3	16,000	4.6 $\times$ 15
				6.0 $\times$ 15
TSKgel SuperHZ2500	$2 \times 10^4$	3	16,000	4.6 $\times$ 15
				6.0 $\times$ 15
TSKgel SuperHZ3000	$6 \times 10^4$	3	16,000	4.6 $\times$ 15
				6.0 $\times$ 15
TSKgel SuperHZ4000	$4 \times 10^5$	3	16,000	4.6 $\times$ 15
				6.0 $\times$ 15
TSKgel SuperH2M-N	$7 \times 10^5$	3	16,000	4.6 $\times$ 15
				6.0 $\times$ 15
TSKgel SuperH2M-M	$4 \times 10^6$	3&5	16,000	4.6 $\times$ 15
				6.0 $\times$ 15
TSKgel SuperH2M-H	$4 \times 10^8$ (Estimate)	10	9,000	4.6 $\times$ 15
				6.0 $\times$ 15

**Table-2**

Feature	Advantage
1) Adoption of fine particles for packing material	<ul style="list-style-type: none"> <li>Short measurement time is achieved. Resolution equivalent to conventional columns (30cm) can be obtained in 1/2 measurement time.</li> <li>Resolution does not deteriorate even under a high flow rate.</li> </ul>
2) Semi-micro column	<ul style="list-style-type: none"> <li>Reduction in solvent consumption (running costs, effluent processing costs) 1/6 to 1/3 solvent consumption compared to conventional columns.</li> </ul>
3) Optimization of particle size in the packing material	<ul style="list-style-type: none"> <li>Shear degradation in polymers with high molecular weight can be prevented.</li> </ul>
4) Adoption of low-absorption packing materials	<ul style="list-style-type: none"> <li>Applicable to wide range of samples.</li> </ul>

## 3. Basic Properties


### 3-1. Resolution

As shown in **Table-1**, columns for separation of low-molecular weight substances such as oligomers, TSKgel SuperHZ1000 - SuperHZ4000 and TSKgel SuperHZM-N that is required high resolution are packed with packing materials which particle size is 3 $\mu$ m, and have high theoretical plates values (per unit length) which are twice of those of the H<sub>XL</sub> series. As a result, equal resolution can be obtained with half the measurement time of the conventional products as shown in **Figure-3**. To sufficiently deliver the performance of such high-performance columns, optimization in the detector, sample injector, tubing in the system, etc. (dead volume reduction) is inevitable in addition to optimization of measurement conditions. Please use due caution when using the semi-micro column with inner diameter of 4.6mm, in which the effect becomes significant (column performance can be sufficiently delivered by using HLC-

8220GPC. In addition, micro flow cell is required for UV detector).


As shown in **Figure-2**, the advantage of TSKgel SuperHZM series (HZ mix series) is that it has a wide separation range of molecular weights with linearity in the curve. That is, it has a wide range of measurable molecular weights with a column of 1 grade. Furthermore, similar to the SuperHZ series for oligomer measurement, TSKgel SuperHZM-N (molecular weight separation range  $7 \times 10^5$  - 266), SuperHZM-M (molecular weight separation range  $4 \times 10^6$  - 266), and SuperHZM-H (molecular weight separation range  $1 \times 10^8$  - 1,000 (estimate)) are capable of measuring molecular weight and molecular weight distribution in half the measurement time of the conventional products.

**Figure-4** shows a comparison of polyisobutylene chromatograms by GMH<sub>XL</sub>, which is a packing material with particle size of 9 $\mu$ m, and SuperHZM-H with particle size of 10 $\mu$ m.

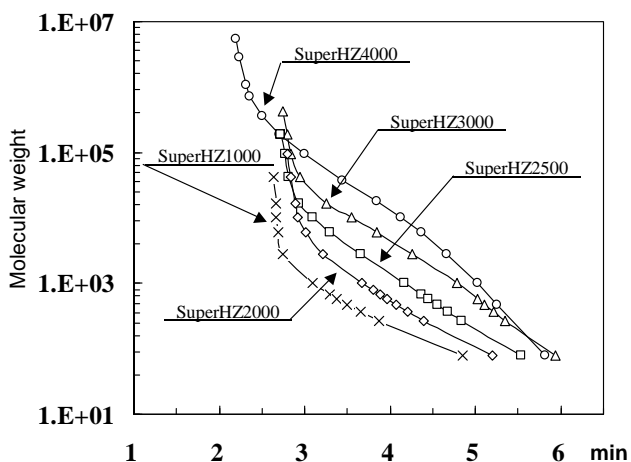

**Figure-1 Calibration Curves of TSKgel SuperHZ Series**

Column: TSKgel SuperHZ series  
(4.6mm I.D.  $\times$  15cm)  
Eluent: THF  
Flow rate: 0.35mL/min  
Temperature: 25°C  
Sample: standard polystyrene  
Injection volume: 2 $\mu$ L

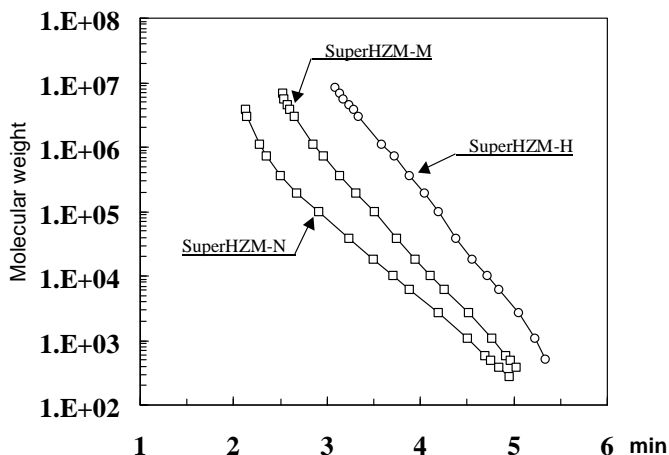


**Figure-2 Calibration Curves of TSKgel SuperHZM Series**

Column: TSKgel SuperHZM series  
(4.6mm I.D.  $\times$  15cm)  
Eluent: THF  
Flow rate: 0.35mL/min  
Temperature: 25°C  
Sample: standard polystyrene  
Injection volume: 2 $\mu$ L

### 3-2. Dependence of Height Equivalent to a Theoretical Plate (HETP) on Flow Rate

**Figure-5** shows the results of comparing the flow rate dependence of HETP in TSKgel G2500H<sub>XL</sub> (particle size 5 $\mu$ m) and TSKgel SuperHZ2500 (particle size 3 $\mu$ m) using a low-molecular-weight compound as sample. While HETP gradually increases in the high flow rate region in G2500H<sub>XL</sub> with larger particle size, little change in HETP is seen in the high flow rate region in SuperHZ2500 with small particle size, indicating that its column efficiency does not deteriorate even under a high flow rate. HETP flow rate dependency of each of SuperHZ2500, SuperH2M-N, SuperH2M-M, and SuperH2M-H (inner diameter 4.6mm and 6.0mm) columns is shown in **Figures-6, -7, -8, and -9**, respectively.


In actual measurement, HETP flow rate dependency will also depend on the sample's molecular size (molecular weight), eluent type (viscosity) and measurement temperature. Especially with high molecular weight polymers, HETP gradually increases as the flow rate increases and this tendency becomes more significant as the molecular weight becomes larger. Therefore, it can be generally said that low flow rate is desired when a high molecular weight polymer is measured.

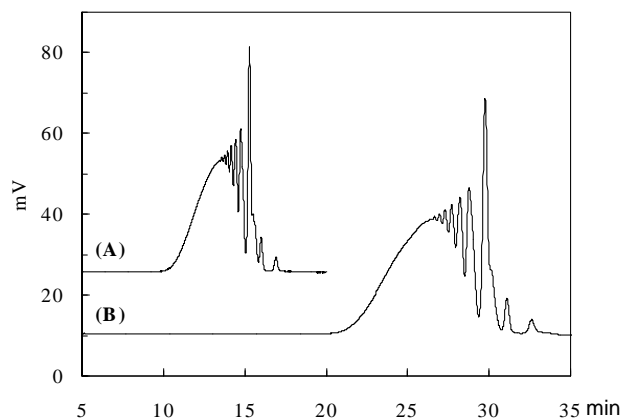

**Figure-3 Comparison between SuperHZ and GMH<sub>XL</sub>**


Column: (A) TSKgel SuperHZ  
(4000 + 3000 + 2500)  
(4.6mm I.D. × 15cm × 3)  
(B) TSKgel Super G  
(4000 + 3000 + 2500) H<sub>XL</sub>  
(7.8mm I.D. × 30cm × 3)

Eluent: THF  
Flow rate: (A) 0.35mL/min  
(B) 1.0mL/min  
Temperature: 40°C  
Detection: RI  
Sample: phenolic resin  
Injection volume: (A) 5 $\mu$ L, (B) 30 $\mu$ L

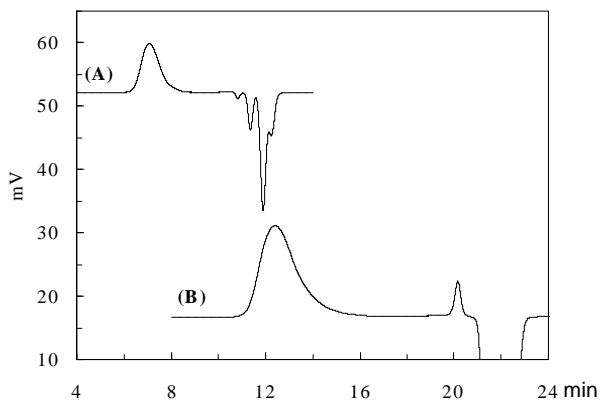


**Figure-4 Comparison between SuperH2M-H and GMH<sub>XL</sub>**

Column: (A) TSKgel SuperH2M-H  
(4.6mm I.D. × 15cm × 2)  
(B) TSKgel GMH<sub>XL</sub>  
(7.8mm I.D. × 30cm × 2)  
Eluent: THF  
Flow rate: (A) 0.35mL/min  
(B) 1.0mL/min  
Temperature: 40°C  
Detection: RI  
Sample: phenolic resin  
Injection volume: (A) 10 $\mu$ L, (B) 100 $\mu$ L

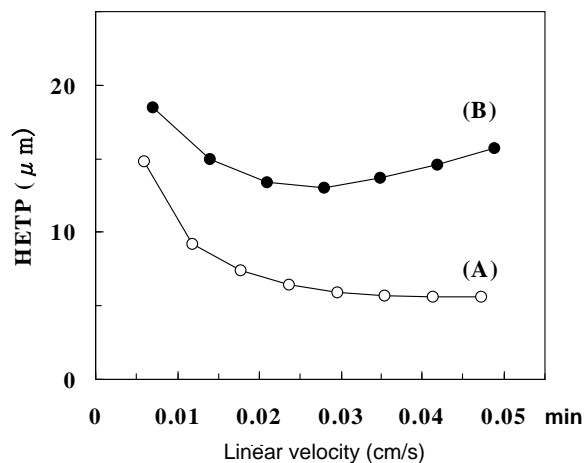


**Figure-5 Comparison of the Effect of Linear Velocity on HETP in SuperHZ2500 and G2500H<sub>XL</sub>**

Column: (A) TSKgel SuperHZ2500  
(6.0mm I.D. × 15cm)  
(B) TSKgel G2500H<sub>XL</sub>  
(7.8mm I.D. × 30cm)  
Eluent: THF  
Temperature: 25°C  
Sample: (A) dicyclohexyl phthalate (DCHP),  
(B) benzene  
Injection volume: (A) 2 $\mu$ L  
(B) 20 $\mu$ L

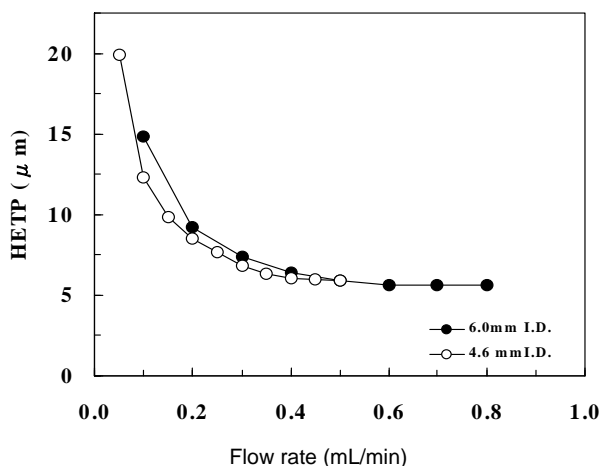


**Figure-6 Relationship between HETP and Flow Rate in SuperHZ2500**

Column: TSKgel SuperHZ2500  
 (6.0mm I.D. × 15cm)  
 (4.6mm I.D. × 15cm)  
 Eluent: THF  
 Temperature: 25°C  
 Sample: DCHP  
 Injection volume: 2μL (6.0mm I.D. × 15cm)  
 1μL (4.6mm I.D. × 15cm)

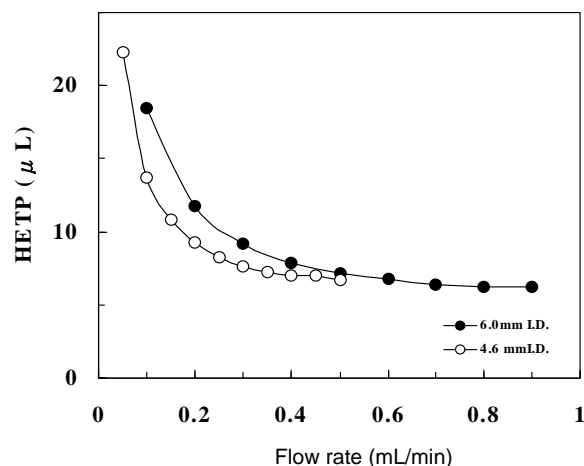


**Figure-7 Relationship between HETP and Flow Rate in SuperH2M-N**

Column: TSKgel SuperH2M-N  
 (6.0mm I.D. × 15cm)  
 (4.6mm I.D. × 15cm)  
 Eluent: THF  
 Temperature: 25°C  
 Sample: DCHP  
 Injection volume: 2μL (6.0mm I.D. × 15cm)  
 1μL (4.6mm I.D. × 15cm)

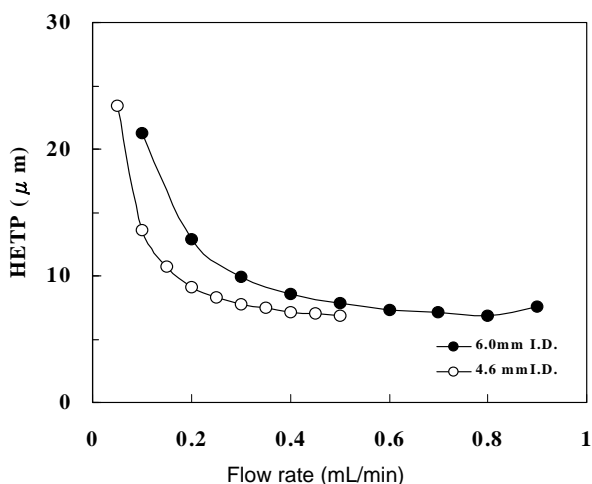

**Figure-8 Relationship between HETP and Flow Rate in SuperH2M-M**

Column: TSKgel SuperH2M-M  
 (6.0mm I.D. × 15cm)  
 (4.6mm I.D. × 15cm)  
 Eluent: THF  
 Temperature: 25°C  
 Sample: DCHP  
 Injection volume: 2μL (6.0mm I.D. × 15cm)  
 1μL (4.6mm I.D. × 15cm)

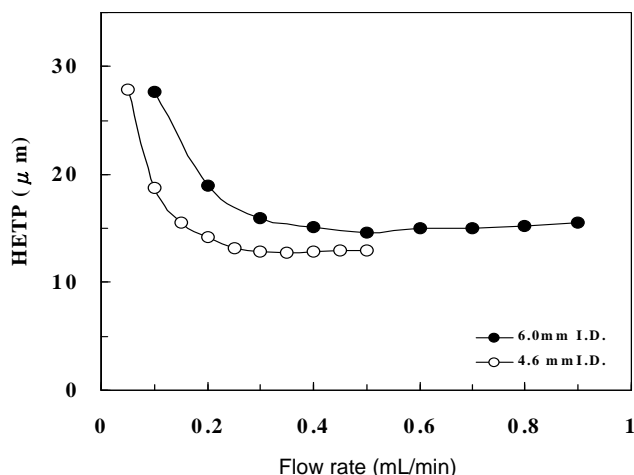

**Figure-9 Relationship between HETP and Flow Rate in SuperH2M-H**


Column: TSKgel SuperH2M-H  
 (6.0mm I.D. × 15cm)  
 (4.6mm I.D. × 15cm)  
 Eluent: THF  
 Temperature: 25°C  
 Sample: DCHP  
 Injection volume: 2μL (6.0mm I.D. × 15cm)  
 1μL (4.6mm I.D. × 15cm)

### 3-3. Effect of Sample Injection Volume

Since TSKGEL SuperHZ series has small packing material particle size as well as small column size, sample injection volume becomes an important factor for obtaining sufficient column performance. As shown in the case of SuperHZ2500 in **Figure-10**, HETP value grows larger and the column performance deteriorates as the sample injection volume increases. This HETP deterioration due to increase in the sample injection volume is seen with more significance in high-performance columns with small packing material particle sizes. As shown in **Figures-10, -11, -12, and -13**, the maximum sample injection volume of each column is 2 $\mu$ L (inner diameter 4.6mm I.D.) and 4 $\mu$ L (inner diameter 6.0mm I.D.) for SuperHZ2500 and SuperHBM-N for oligomer measurement, 2 $\mu$ L (inner diameter 4.6mm I.D.)


and 5 $\mu$ L (inner diameter 6.0mm I.D.) for Super HZM-M for synthetic polymer measurement, and 5 $\mu$ L (inner diameter 4.6mm I.D.) and 10 $\mu$ L (inner diameter 6.0mm I.D.) for SuperHBM-H.

In measurement of polymers including oligomers in which high resolution is most important, it is necessary that measurement is taken with high sample concentration so that the sample injection volume is less than the maximum volume. On the other hand, for synthetic polymers that do not contain oligomers and have relatively large molecular weights, measurement should be taken with low concentration and large injection volume because the dependency of measured molecular weight on the sample injection volume is small. At this time, sample injection volume should be optimized with consideration of overloading, detection sensitivity, etc.

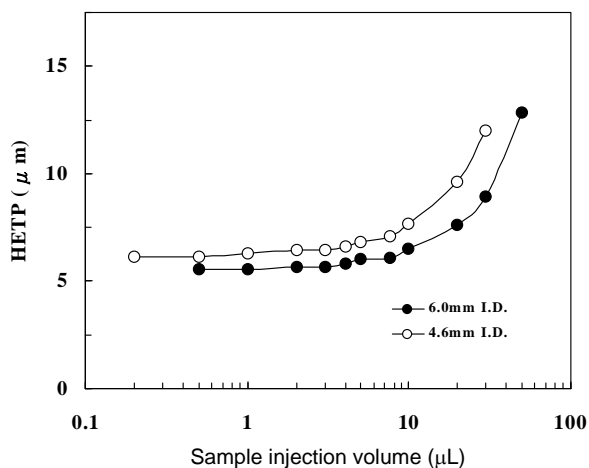


**Figure-10 Relationship between HETP and Sample Injection Volume in SuperHZ2500**

Column: TSKgel SuperHZ2500  
 (6.0mm I.D. × 15cm)  
 (4.6mm I.D. × 15cm)  
 Eluent: THF  
 Flow rate: 0.6mL/min (6.0mm I.D. × 15cm)  
 0.35mL/min (4.6mm I.D. × 15cm)  
 Temperature: 25°C  
 Sample: DCHP

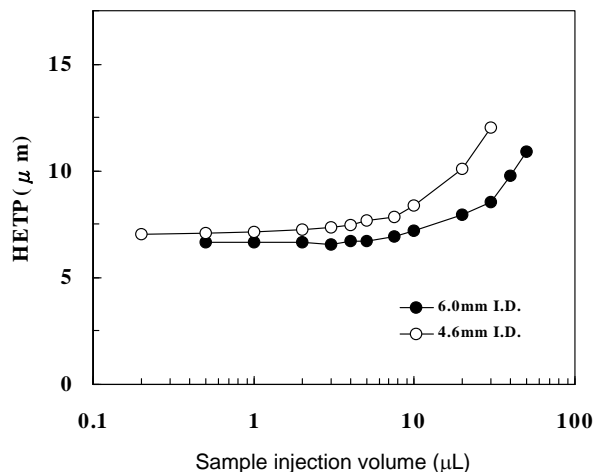


**Figure-11 Relationship between HETP and Sample Injection Volume in SuperHBM-N**

Column: TSKgel SuperHBM-N  
 (6.0mm I.D. × 15cm)  
 (4.6mm I.D. × 15cm)  
 Eluent: THF  
 Flow rate: 0.6mL/min (6.0mm I.D. × 15cm)  
 0.35mL/min (4.6mm I.D. × 15cm)  
 Temperature: 25°C  
 Sample: DCHP

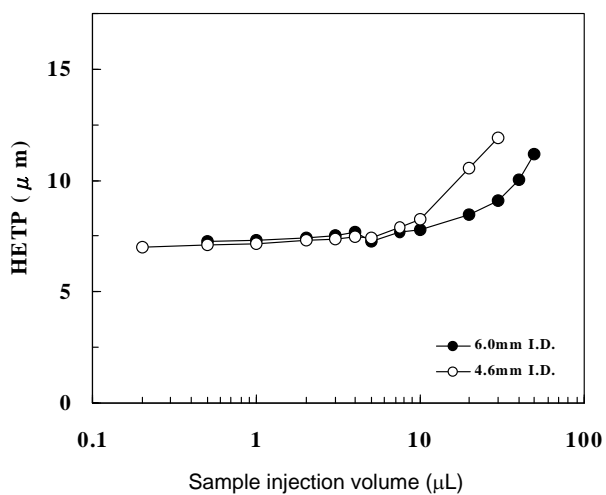

**Figure-12 Relationship between HETP and Sample Injection Volume in SuperHZM-M**

Column: TSKgel SuperHZM-M  
 (6.0mm I.D. × 15cm)  
 (4.6mm I.D. × 15cm)  
 Eluent: THF  
 Flow rate: 0.6mL/min (6.0mm I.D. × 15cm)  
 0.35mL/min (4.6mm I.D. × 15cm)  
 Temperature: 25°C  
 Sample: DCHP

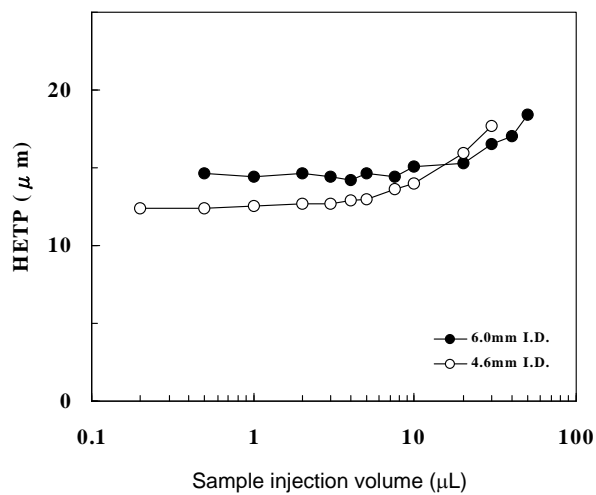

**Figure-13 Relationship between HETP and Sample Injection Volume in SuperHZM-H**

Column: TSKgel SuperHZM-H  
 (6.0mm I.D. × 15cm)  
 (4.6mm I.D. × 15cm)  
 Eluent: THF  
 Flow rate: 0.6mL/min (6.0mm I.D. × 15cm)  
 0.35mL/min (4.6mm I.D. × 15cm)  
 Temperature: 25°C  
 Sample: DCHP


### 3-4. Effect of Sample Concentration

Since the apparent sample molecule size becomes small when the sample concentration increases, the overload phenomenon in which elution volume becomes slow starts to appear. As shown in **Figure-14**, the overload phenomenon is similar to the effect of sample injection volume. It becomes more significant as the particle size of the packing material becomes smaller and slowing in elution volume is observed from lower concentrations as the molecular weight of the measured sample is larger. Caution is required since the calibration curve will appear as shown in **Figure-15** and the molecular weight of the measured sample will be calculated larger when calibration curve is obtained under such circumstances. Proper molecular weight measurement will be difficult when overload phenomenon occurs in either of the molecular weight standard for calibration curve or in the measured sample.

The relationship between the standard polystyrene sample concentration and elution time in TSKGEL SuperHZ mix grade (4.6mm I.D.) is shown in **Figures-16**, **-17** and **-18**. When actual polymers are measured, optimal sample concentration varies depending on the molecular weight distribution and the column to be used in the molecular weight measurement. However, when epoxy


resin is measured using SuperH2M-N (4.6mm I.D.) as shown in **Figure-19**, little change is seen in various average molecular weights and resolutions as long as the sample concentration stays below 20g/L (load 100 $\mu$ g), which is relatively high, since the molecular weight of the measured sample is small and the effect of overloading is not seen. **Figure-20** shows the chromatograms of epoxy resin under different sample concentrations. With concentration of 100g/L (load 500 $\mu$ g), slowing in elution position due to overloading is seen.

**Figures-21**, and **-22** show the relationship between the sample concentration and molecular weight. Measurement needs to be performed with load of 20 $\mu$ g or less when polystyrene SRM706 (Mw257000) is used as the measured sample in SuperH2M-M (4.6mm I.D.), because deterioration in molecular weight value due to overloading is seen at sample concentration of 2g/L (load 20 $\mu$ g) or higher. On the other hand, when polyisobutylene is measured in SuperH2M-H (4.6mm I.D.), deterioration in the molecular weight value is seen from the sample concentration of 1g/L (load 10 $\mu$ g) because the molecular weight of the sample is large. Thus it is necessary that measurement must be taken under a sufficiently low concentration when a sample with large molecular weight is measured.

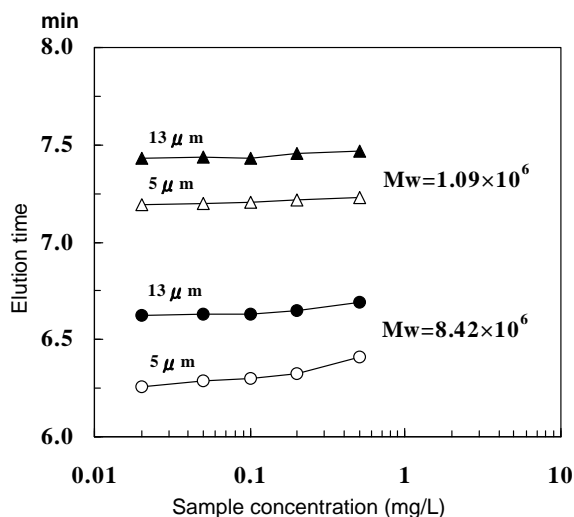

**Figure-14 Effect of Sample Concentration on Elution Time**



Column: TSKgel SuperH2M-H × 2  
(4.6mm I.D. × 15cm × 2)  
Eluent: THF  
Flow rate: 0.35mL/min  
Temperature: 40°C  
Detection: RI  
Sample: standard polystyrene  
Injection volume: 10 $\mu$ L

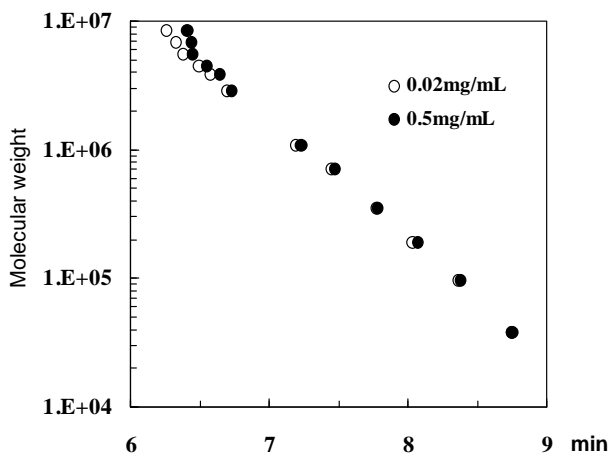

**Figure-15 Effect of Sample Concentration on Elution Time (Calibration Curve)**



Column: TSKgel SuperH2M-H × 2  
(4.6mm I.D. × 15cm × 2)  
Eluent: THF  
Flow rate: 0.35mL/min  
Temperature: 40°C  
Detection: RI  
Sample: standard polystyrene  
Injection volume: 10 $\mu$ L



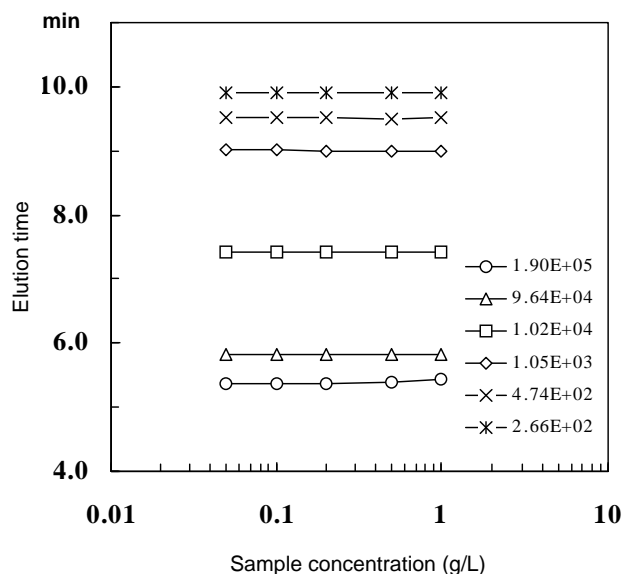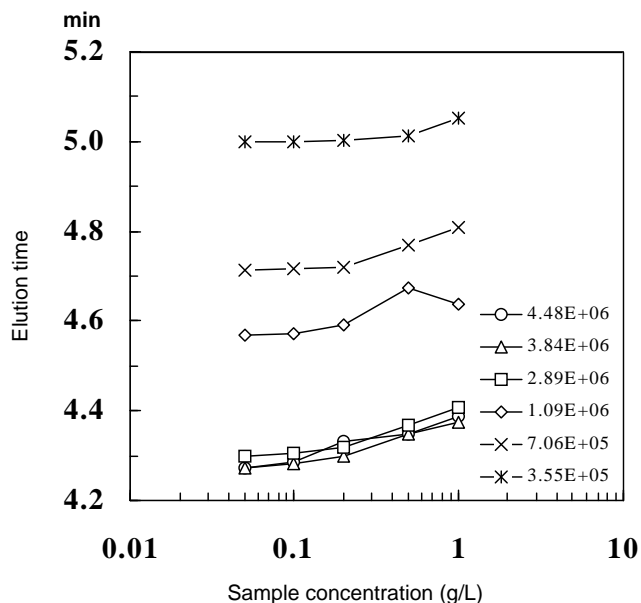



**Figure-16 Effect of Sample Concentration on Elution Time**

Column: TSKgel SuperH2M-N × 2  
(4.6mm I.D. × 15cm × 2)  
Eluent: THF  
Flow rate: 0.35mL/min  
Temperature: 40°C  
Detection: RI  
Sample: standard polystyrene  
Injection volume: 5µL

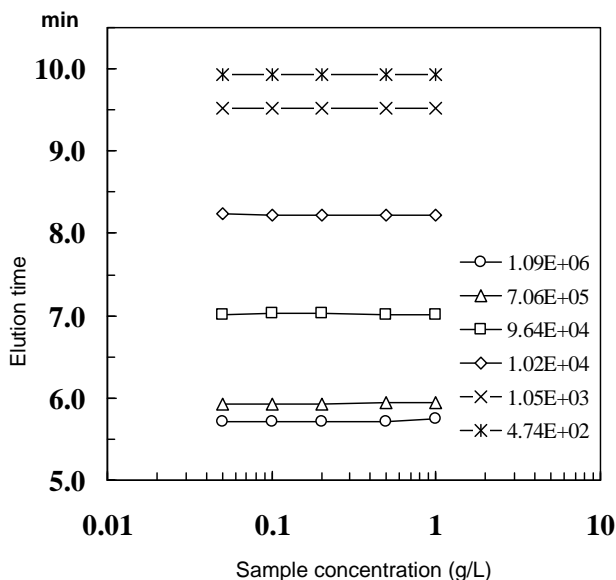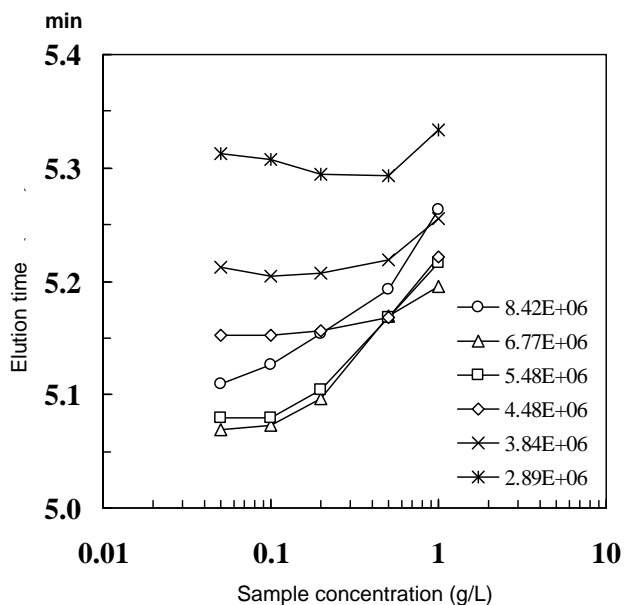


**Figure-17 Effect of Sample Concentration on Elution Time**

Column: TSKgel SuperH2M-M × 2  
(4.6mm I.D. × 15cm × 2)  
Eluent: THF  
Flow rate: 0.35mL/min  
Temperature: 40°C  
Detection: RI  
Sample: standard polystyrene  
Injection volume: 10µL

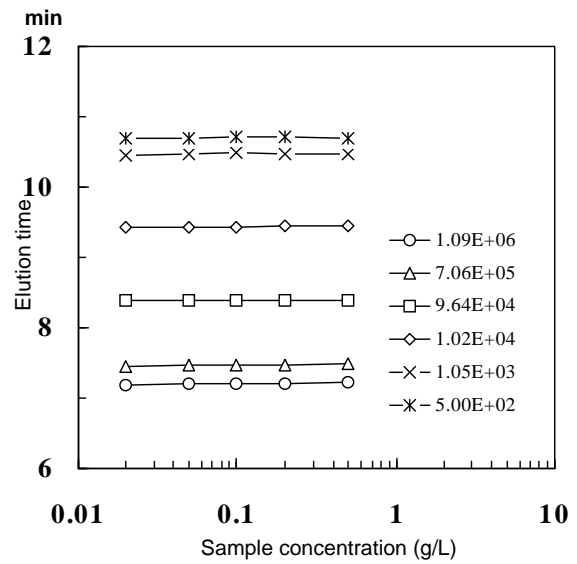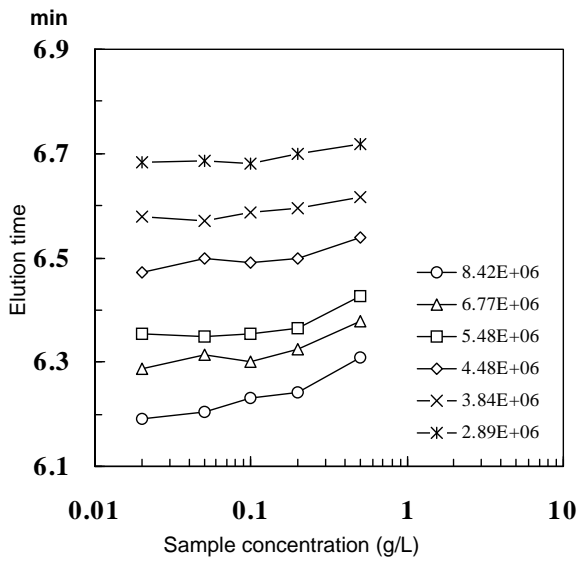


**Figure-18 Effect of Sample Concentration on Elution Time**

Column: TSKgel SuperHZM-H × 2  
(4.6mm I.D. × 15cm × 2)  
Eluent: THF  
Flow rate: 0.35mL/min  
Temperature: 40°C  
Detection: RI  
Sample: standard polystyrene  
Injection volume: 5μL

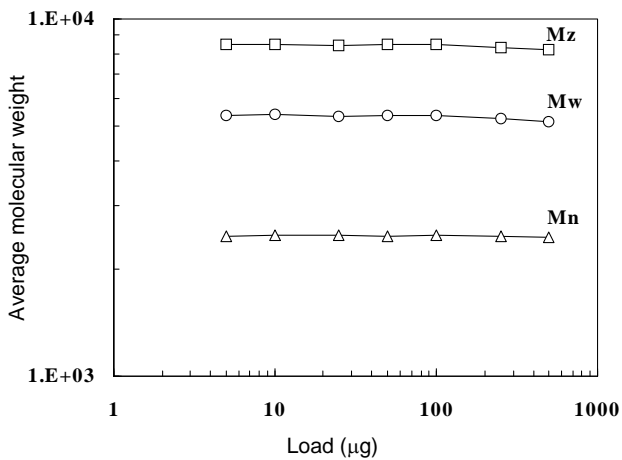


**Figure-19 Effect of Sample Injection Volume on Average Molecular Weight**

Column: TSKgel SuperHZM-N  
(4.6mm I.D. × 15cm × 2)  
Eluent: THF  
Flow rate: 0.35mL/min  
Temperature: 40°C  
Detection: RI  
Sample: epoxy resin  
Injection volume: 5μL

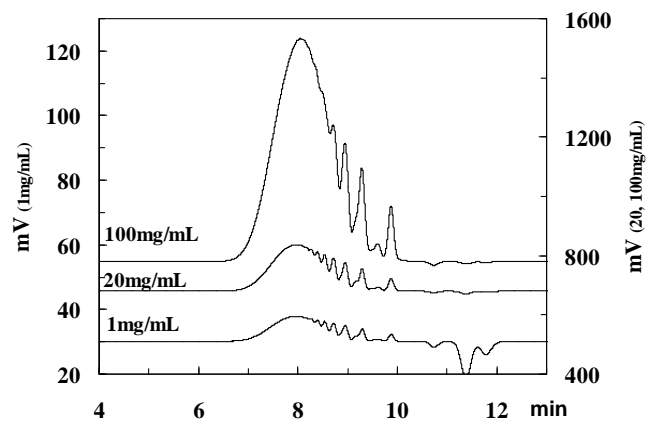


**Figure-20 Effect of Sample Load on Chromatogram**

Column: TSKgel SuperHZM-N  
(4.6mm I.D. × 15cm × 2)  
Eluent: THF  
Flow rate: 0.35mL/min  
Temperature: 40°C  
Detection: RI  
Sample: epoxy resin  
Injection volume: 5μL

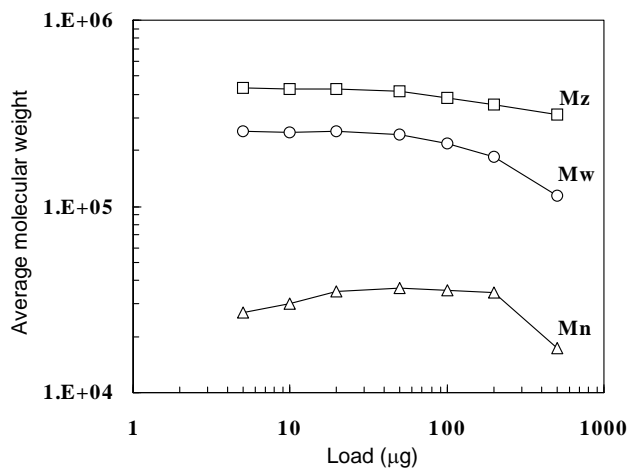

**Figure-21 Effect of Sample Injection Volume on Average Molecular Weight**

Column: TSKgel SuperH2M-M  
 (4.6mm I.D. × 15cm × 2)  
 Eluent: THF  
 Flow rate: 0.35mL/min  
 Temperature: 40°C  
 Detection: RI  
 Sample: polystyrene SRM706  
 Injection volume: 10µL

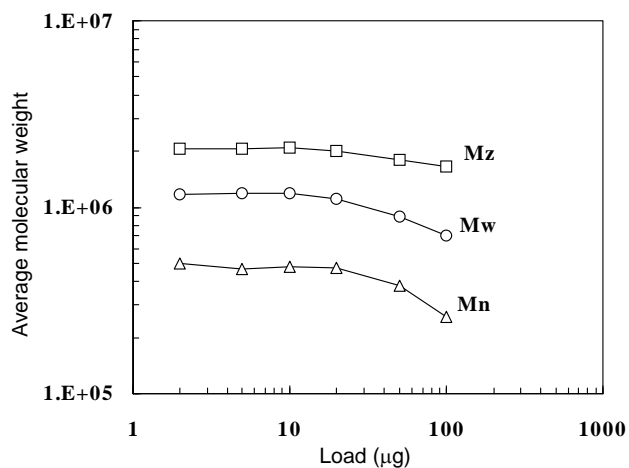

**Figure-22 Effect of Sample Injection Volume on Average Molecular Weight**


Column: TSKgel SuperH2M-H  
 (4.6mm I.D. × 15cm × 2)  
 Eluent: THF  
 Flow rate: 0.35mL/min  
 Temperature: 40°C  
 Detection: RI  
 Sample: polyisobutylene  
 Injection volume: 10µL

### 3-5. Optimization of Packing Material Particle Size

When measurement is taken using a column packed with fine-particle packing materials or under a high flow rate, it becomes more likely that chain breaking in polymers occur. Therefore, in TSKGEL SuperHZ series, optimal particle size that corresponds to the range of molecular weight fractionation is adopted in each of TSKgel SuperH2M-M and SuperH2M-H, which are the columns for polymers.

**Figure-23** shows the calibration curve obtained by measuring the standard polystyrene using packing materials of different particle sizes and different measurement flow rates. When a packing material with particle size of 3 $\mu$ m is used for samples with molecular weight of 1 million or higher, delay in elution due to breaking in branched chains is seen at 0.35mL/min. Furthermore, delay in elution can also be seen for samples with molecular weight of several hundred thousands or larger. This phenomenon appears more

drastically as the particle size becomes smaller and the measurement flow rate becomes higher. On the other hand, when a packing material with particle size of 13 $\mu$ m is used, delay in sample elution is not seen throughout the molecular weight range under any flow rate. Based on this fact, SuperH2M-H, which targets measurement of average molecular weight of several millions to several hundred thousands, has been packed with a packing material of particle size 10 $\mu$ m and SuperH2M-M, which targets measurement of average molecular weight of several hundred thousands to several ten thousands, has been packed with a packing material of particle size 5 or 3 $\mu$ m. In addition, a packing material with particle size of 3 $\mu$ m has been used for SuperH2M-N, in which high resolution is focused to target samples that include oligomers and have average molecular weight of several ten thousands. As you can see, SuperH2M series has been packed with packing materials with optimal particle sizes that correspond to the molecular weights of the target samples.

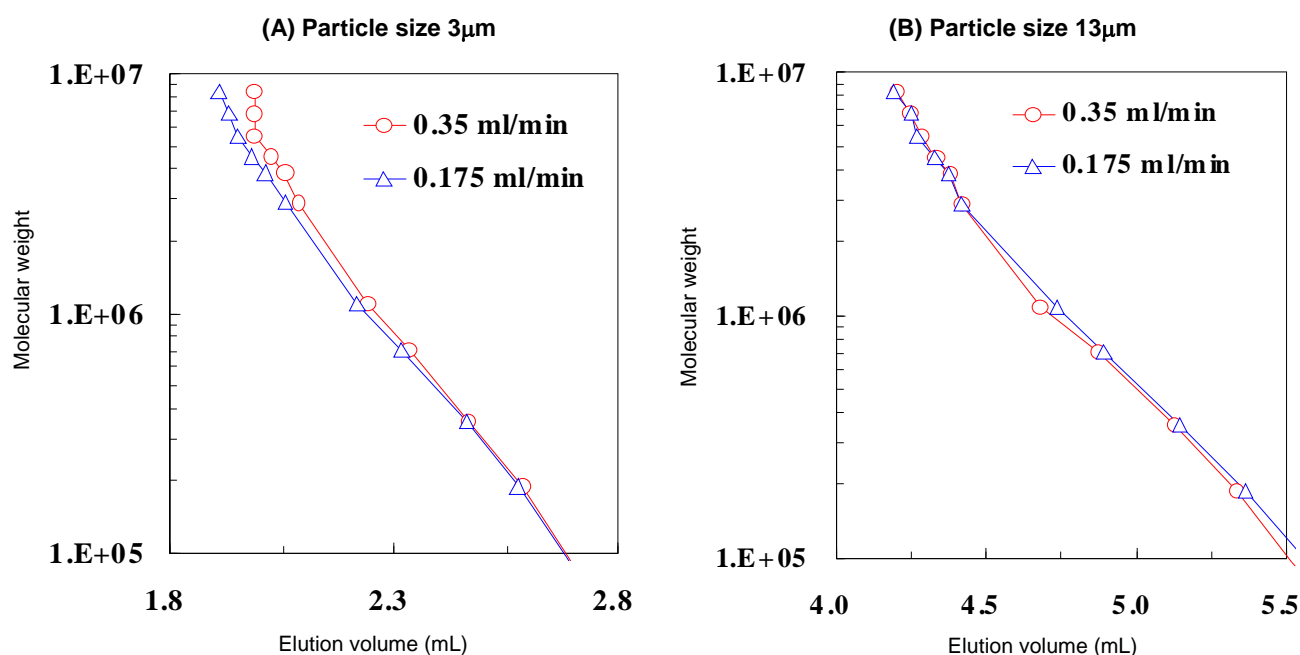


**Figure-23 Effect of Measurement Flow Rate and Particle Size on Calibration Curve**

Column: styrene-divinylbenzene copolymerization gel  
 (A): (4.6mm I.D.  $\times$  15cm)  
 (B): (4.6mm I.D.  $\times$  30cm)  
 Eluent: THF  
 Flow rate: 0.35mL/min  
 Temperature: 40°C  
 Detection: RI  
 Sample: standard polystyrene  
 Injection volume: (A) 5 $\mu$ L, (B)10 $\mu$ L


## 4. Application

Figures-24 to -35 show examples of analysis using various samples.

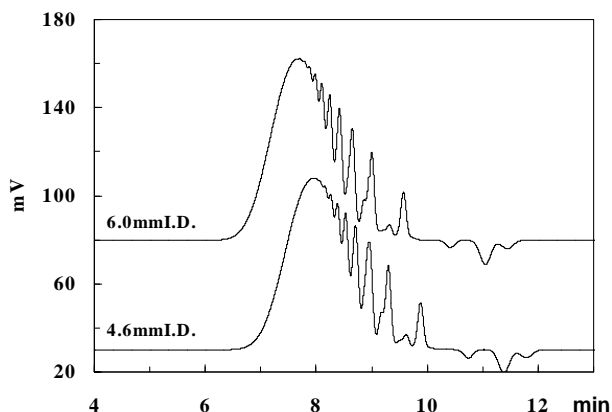


**Figure-24 Chromatogram of Epoxy Resin**

Column: TSKgel SuperHZM-N × 2  
 Eluent: THF  
 Flow rate: 0.35mL/min (4.6mm I.D.)  
 0.6mL/min (6.0mm I.D.)  
 Temperature: 40°C Detection: RI  
 Sample: epoxy resin (10g/L)  
 Injection volume: 5µL (4.6mm I.D.)  
 9µL (6.0mm I.D.)

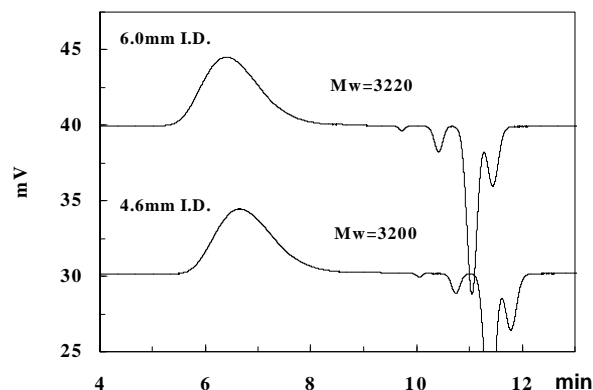


**Figure-25 Chromatogram of Polymethyl Methacrylate**

Column: TSKgel SuperHZM-N × 2  
 Eluent: THF  
 Flow rate: 0.35mL/min (4.6mm I.D.)  
 0.6mL/min (6mm I.D.)  
 Temperature: 40°C Detection: RI  
 Sample: polymethyl methacrylate (1g/L)  
 Injection volume: 5µL (4.6mm I.D.)  
 9µ (6.0mm I.D.)

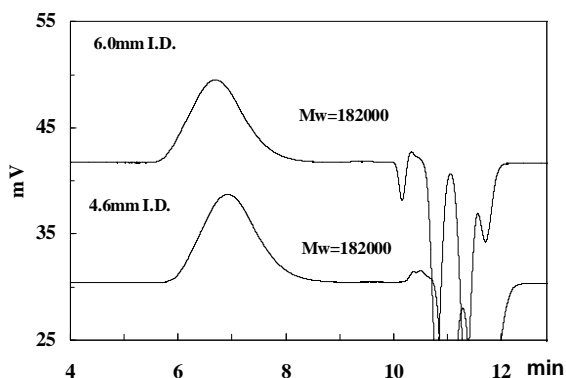


**Figure-26 Chromatogram of Butylmethacrylate-Isobutylmethacrylate Copolymer**

Column: TSKgel SuperHZM-M × 2  
 Eluent: THF  
 Flow rate: 0.35mL/min (4.6mm I.D.)  
 0.6mL/min (6mm I.D.)  
 Temperature: 40°C Detection: RI  
 Sample: butylmethacrylate-isobutylmethacrylate copolymer (1g/L)  
 Injection volume: 10µL (4.6mm I.D.)  
 17µL (6.0mm I.D.)

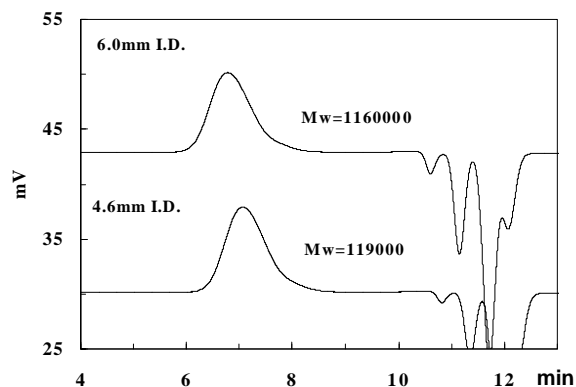


**Figure-27 Chromatogram of Polyisobutylene**


Column: TSKgel SuperHZM-M × 2  
 Eluent: THF  
 Flow rate: 0.35mL/min (4.6mm I.D.)  
 0.6mL/min (6.0mm I.D.)  
 Temperature: 40°C Detection: RI  
 Sample: polyisobutylene (0.5g/L)  
 Injection volume: 10µL (4.6mm I.D.)  
 17µL (6.0mm I.D.)

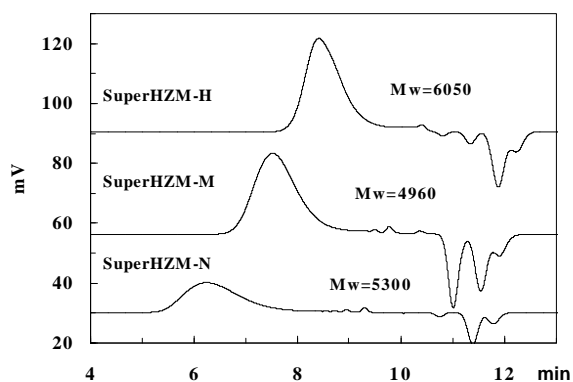


**Figure-28 Chromatogram of Polysulfone**  
 Column: TSKgel SuperHZM × 2 (4.6mm I.D.)  
 Eluent: THF  
 Flow rate: 0.35mL/min  
 Temperature: 40°C Detection: RI  
 Sample: polysulfone (1g/L)  
 Injection volume: 5µL (SuperHZM-N),  
 10µL (SuperHZM-M, -H)

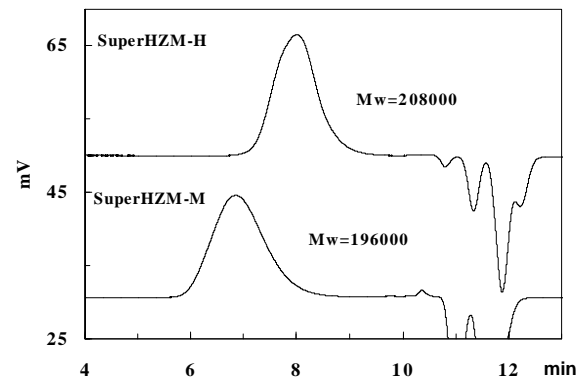


**Figure-29 Chromatogram of 1, 2-Polybutadiene**  
 Column: TSKgel SuperHZM × 2 (4.6mm I.D. × 15cm × 2)  
 Eluent: THF  
 Flow rate: 0.35mL/min  
 Temperature: 40°C Detection: RI  
 Sample: polybutadiene (1.0g/L)  
 Injection volume: 10µL

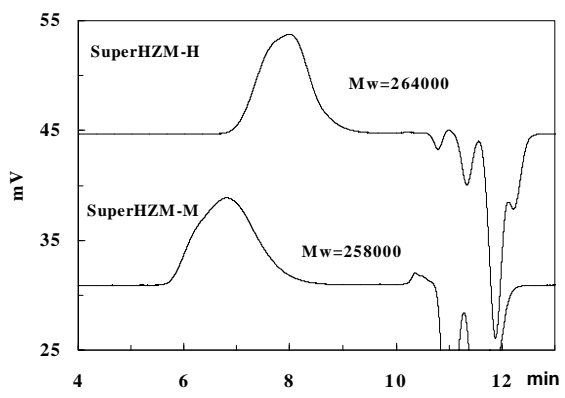

**Figure-30 Chromatogram of Polyethylmethacrylate**  
 Column: TSKgel SuperHZM × 2 (4.6mm I.D. × 15cm × 2)  
 Eluent: THF  
 Flow rate: 0.35mL/min  
 Temperature: 40°C Detection: RI  
 Sample: polyethylmethacrylate (1.0g/L)  
 Injection volume: 10µL

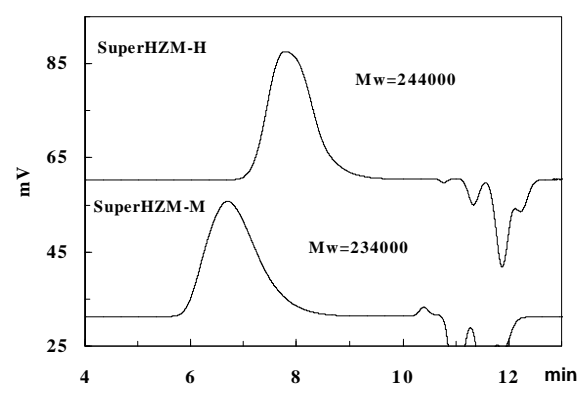


**Figure-31 Chromatogram of Polystyrene**  
 Column: TSKgel SuperHZM × 2 (4.6mm I.D. × 15cm × 2)  
 Eluent: THF  
 Flow rate: 0.35mL/min  
 Temperature: 40°C Detection: RI  
 Sample: polystyrene (1.0g/L)  
 Injection volume: 10µL

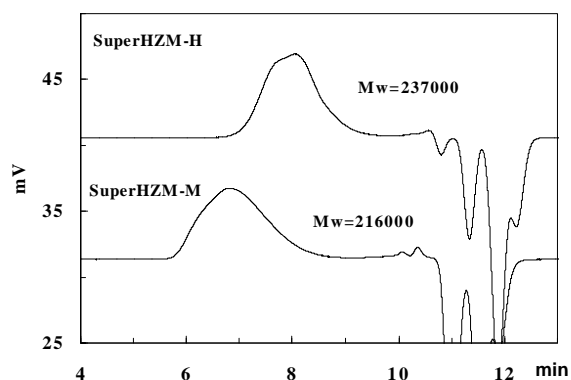


**Figure-32 Chromatogram of Polyvinyl Acetate**

Column: TSKgel SuperH2M × 2  
(4.6mm I.D. × 15cm × 2)  
Eluent: THF  
Flow rate: 0.35mL/min  
Temperature: 40°C Detection: RI  
Sample: polyvinyl acetate (1.0g/L)  
Injection volume: 10µL

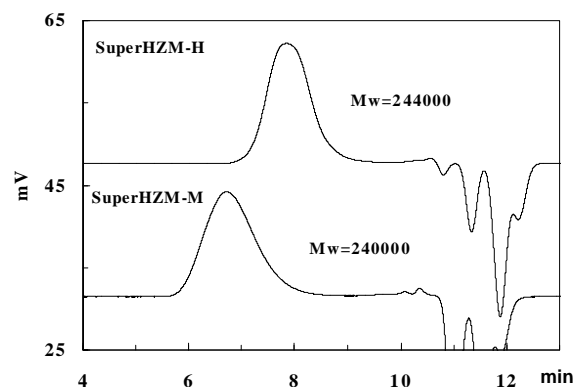


**Figure-33 Chromatogram of Polyvinyl Chloride**

Column: TSKgel SuperH2M × 2  
(4.6mm I.D. × 15cm × 2)  
Eluent: THF  
Flow rate: 0.35mL/min  
Temperature: 40°C Detection: RI  
Sample: polyvinyl chloride (1.0g/L)  
Injection volume: 10µL

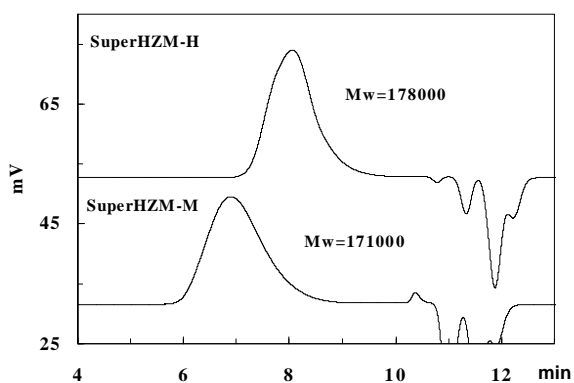

**Figure-34 Chromatogram of Styrene-Acrylonitrile Copolymer**

Column: TSKgel SuperH2M × 2  
(4.6mm I.D. × 15cm × 2)  
Eluent: THF  
Flow rate: 0.35mL/min  
Temperature: 40°C Detection: RI  
Sample: styrene-acrylonitrile copolymer (1.0g/L)  
Injection volume: 10µL

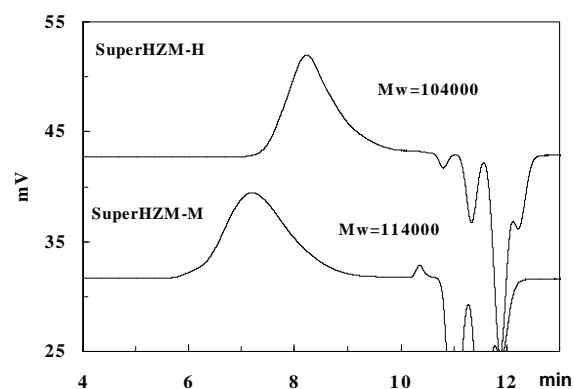

**Figure-35 Chromatogram of Vinyl Alcohol-Vinyl Butyral Copolymer**

Column: TSKgel SuperH2M × 2  
(4.6mm I.D. × 15cm × 2)  
Eluent: THF  
Flow rate: 0.35mL/min  
Temperature: 40°C Detection: RI  
Sample: vinyl alcohol-vinyl butyral copolymer (1.0g/L)  
Injection volume: 10µL

## 5. Conclusion

By optimizing the packing material particle size for each molecular weight measurement range, TSKGEL SuperHZ series is capable of performing measurement under the analysis time that is half of the conventional H<sub>XL</sub> series while maintaining the high resolution in low-molecular weight range and for oligomers and eliminating the need to worry about chain breaking in the high-molecular weight range.

As described earlier in the main text, it is important that measurement is taken by an optimized system under the optimal measurement conditions in order to sufficiently deliver the performance of a high-performance column. While columns with inner diameters 4.6mm and 6.0mm are available, we recommend that measurement be taken in the high-speed GPC system, HLC-8220GPC, which is not influenced by 25C change, etc. and excels in liquid feeding repeatability, when columns with inner diameter of 4.6mm for low flow rate are used.