



TOSOH

No.106

SEPARATION REPORT

Aqueous SEC Columns for Analysis of Cationic Polymers: TSKgel PW_{XL}-CP Series

— Contents —

	Page
1. Introduction	1
2. Features	1
3. Basic characteristics	2
3-1. Pore characteristics	2
3-2. Sample injection volume and height equivalent to theoretical plate (HETP)	2
3-3. Dependence of height equivalent to theoretical plate (HETP) on flow rate	3
3-4. Separation performance	3
4. Dependence of calibration curves on eluents	4
4-1. Temperature dependence	4
4-2. Dependence on pH in eluent	4
4-3. Effects of organic solvent concentration in solvent	5
4-4. Flow rate dependency	5
5. Elution characteristics of cationic polymers	6
5-1. Comparison of elution using commercial SEC column	6
5-2. Measurement repeatability	8
5-3. Application examples analyzing various cationic polymers	9
6. Analysis of absolute molecular mass of cationic polymers by SEC-MALS	10
6-1. Absolute molecular mass of polydiallyl dimethyl ammonium chloride	10
6-2. Absolute molecular mass of PDADM-NH ₄ Cl copolymer	11
7. Conclusion	11

1. Introduction

Polymer and silica matrices are widely used as packing materials in aqueous size exclusion chromatography (SEC) for analyzing molecular mass distributions of water-soluble polymers. At Tosoh, we have developed and marketed the TSKgel PW and PW_{XL} series of aqueous SEC columns in which hydrophilic vinyl polymers are used as matrices.

When analyzing the molecular mass distribution of water-soluble polymers by SEC, it is generally necessary to optimize the type and concentration of the salt, as well as the pH of the eluent, in order to inhibit interaction between the sample and the packing material.

In particular, for SEC analyses of cationic polymers, to inhibit electrostatic interactions with the packing material, it is necessary to use low-pH eluents or eluents with high salt concentrations. However, even under these conditions, electrostatic interactions with the packing material may not be inhibited, and due to the resulting low sample recovery rates and poor repeatability, adequate molecular mass data is not always obtained.

To overcome these problems, we have developed and started to market the TSKgel PW_{XL}-CP series of aqueous SEC columns to enable analysis of cationic polymers with high recovery rates and good repeatability, even under mild conditions in which neutral eluents are used.

In this report the basic characteristics and application examples of the new TSKgel PW_{XL}-CP series of SEC columns will be discussed.

2. Features

When the TSKgel PW_{XL} series is used to analyze cationic polymers by SEC, adsorption may be observed, caused by electrostatic interactions between the sample and carboxyl groups present in trace amounts on the surface of the packing material. This can make it impossible for adequate molecular mass distribution data to be obtained.

The TSKgel PW_{XL}-CP series, developed using a new method of synthesis, uses the same matrix as the TSKgel PW_{XL} series. However, by improving the ionicity on the surface of the packing material, while maintaining the basic characteristics of the existing series of columns, this series provides excellent repeatability and high recovery rates for cationic polymers, even under conditions in which neutral eluents with low salt concentrations are used.

The basic characteristics of the TSKgel PW_{XL}-CP series are shown in **Tables 1 and 2**.

Table 1 Properties of the TSKgel PW_{XL}-CP series

	TSKgel G3000PW _{XL} -CP	TSKgel G5000PW _{XL} -CP	TSKgel G6000PW _{XL} -CP
Packing material matrix	Hydrophilic vinyl polymer	Hydrophilic vinyl polymer	Hydrophilic vinyl polymer
Particle size	7 μm	10 μm	13 μm
Molecular mass exclusion limit (PEO)	100,000	1,000,000	20,000,000*
Molecular mass separation range (PEO and PEG)	200~50,000	400~500,000	1,000~10,000,000

*Estimated

Table 2 Specifications of TSKgel PW_{XL}-CP series

Product name	Number of theoretical plates (TP/column)	Asymmetry factor	Column size (mm I.D. x cm)
TSKgel G3000PW _{XL} -CP	16,000/30 cm	0.7 ~ 1.6	7.8 x 30
TSKgel G5000PW _{XL} -CP	10,000/30 cm	0.7 ~ 1.6	7.8 x 30
TSKgel G6000PW _{XL} -CP	7,000/30 cm	0.7 ~ 1.6	7.8 x 30

Analysis conditions

Eluent: 0.1 mol/L aqueous solution of sodium nitrate
 Detection: RI
 Flow rate: 1.0 mL/min (7.8 mm I.D. x 30 cm)

Temperature: 25 °C
 Sample: Ethylene glycol (5 g/L)
 Injection volume: 20 μL

3. Basic characteristics

3-1. Pore characteristics

The TSKgel PW_{XL}-CP series, as shown in **Tables 1** and **2**, is provided in three grades with different molecular mass separation ranges capable of accommodating low- to high-molecular mass samples.

Figure 1 shows calibration curves produced with standard polyethylene oxide and polyethylene glycol in a 0.1 mL/L aqueous solution of sodium nitrate.

The TSKgel G3000PW_{XL}-CP, the low-molecular mass grade column, is capable of analysis within a molecular mass separation range of 50,000 to 200. A molecular mass separation range of 500,000 to 400 was established for the TSKgel G5000PW_{XL}-CP, and 10,000,000 to 1,000 for the TSKgel G6000PW_{XL}-CP series, the high-molecular mass grade column. As a result, by connecting these columns together according to the molecular mass and molecular mass distribution of the sample, this series can accommodate a wide range of samples from high-molecular mass polymers to low-molecular mass oligomers.

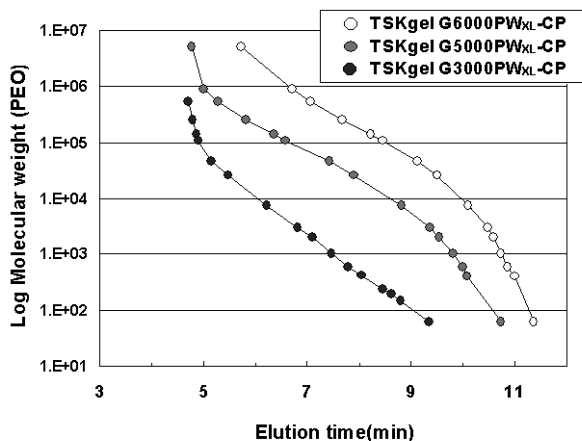


Fig. 1 TSKgel PW_{XL}-CP series calibration curves

Column: TSKgel PW_{XL}-CP series
(7.8 mm I.D. x 30 cm)

Eluent: 0.1 mol/L sodium nitrate

Flow rate: 1.0 mL/min

Detection: RI

Temperature: 25 °C

Samples: Standard polyethylene oxide, polyethylene glycol, and ethylene glycol

Injection volume: 100 μL

3-2. Sample injection volume and height equivalent to theoretical plate (HETP)

Figure 2 shows the dependence of the height equivalent to theoretical plate (HETP) of ethylene glycol on injection volume in the TSKgel PW_{XL}-CP series.

It is clear from this graph that at each grade, if the injection volume exceeds 50 μL when 2 columns are used, HETP begins to increase, and that in the TSKgel PW_{XL}-CP series, the maximum sample injection volume is about 25 μL per column.

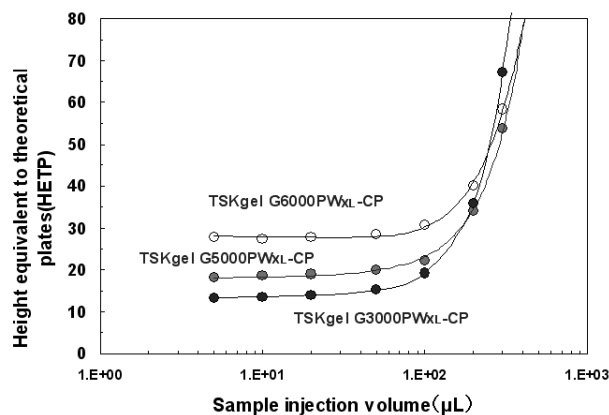


Fig. 2 Relationship between HETP and sample injection volume in TSKgel PW_{XL}-CP series

Column: TSKgel PW_{XL}-CP series
(7.8 mm I.D. x 30 cm) x 2

Eluent: H₂O

Flow rate: 1.0 mL/min

Detection: RI

Temperature: 25 °C

Samples: Ethylene glycol (5 g/L)

Injection volume: 5 to 500 μL

3-3. Dependence of height equivalent to theoretical plate (HETP) on flow rate

Figure 3 confirms the dependence of HETP on flow rate in the TSKgel PW_{XL}-CP series when EG, a low-molecular mass monodisperse material, is used as the sample.

The optimal flow rate (minimum HETP) varies in accordance with the particle size. The optimum flow rate (minimum HETP) with the TSKgel G3000PW_{XL}-CP, which has a small particle size, is relatively high (0.6 to 1.0 mL/min). However, as the particle size increases, the optimal flow rate decreases to 0.5 to 0.8 mL/min with the TSKgel G5000PW_{XL}-CP, and to 0.4 to 0.7 mL/min with the TSKgel G6000PW_{XL}-CP.

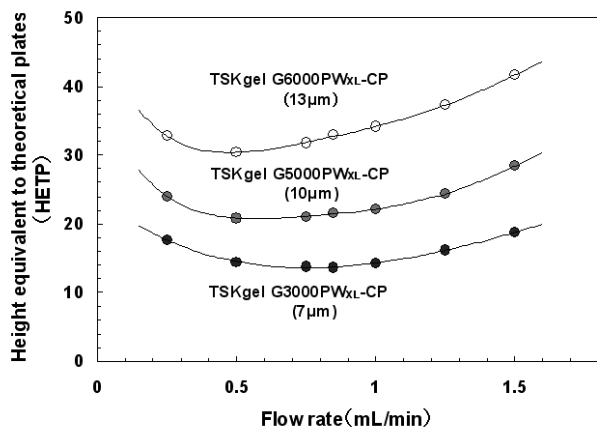


Fig. 3 Relationship between HETP and flow rate in the TSKgel PW_{XL}-CP series

Column: TSKgel PW_{XL}-CP series (7.8 mm I.D. x 30 cm)
 Eluent: H₂O
 Flow rate: 0.25 to 1.5 mL/min
 Detection: RI
 Temperature: 25 °C
 Samples: Ethylene glycol (5 g/L)
 Injection volume: 20 µL

3-4. Separation performance

Although it was previously stated that with regard to pore characteristics, there are no major differences between the TSKgel PW_{XL}-CP series and TSKgel PW_{XL} series of columns, the pore characteristics of the TSKgel G3000PW_{XL}-CP, the low-molecular mass grade column, have been improved in comparison to the TSKgel G3000PW_{XL}.

Fig. 4 shows chromatograms of polyethylene glycol 200 (PEG 200) produced using the TSKgel G3000PW_{XL}-CP and the TSKgel G3000PW_{XL} columns.

The TSKgel G3000PW_{XL}-CP, which has improved pore characteristics in the low-molecular mass range, shows better separation than the TSKgel G3000PW_{XL}.

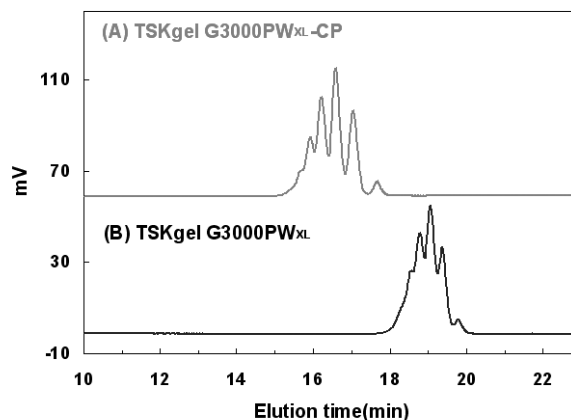


Fig. 4 Separation of Polyethylene Glycol 200 using the TSKgel G3000PW_{XL}-CP and TSKgel G3000PW_{XL}

Column: (A) TSKgel G3000PW_{XL}-CP (7.8 mm I.D. x 30 cm x 2)
 (B) TSKgel G3000PW_{XL} (7.8 mm I.D. x 30 cm x 2)

Eluent: H₂O
 Flow rate: 1.0 mL/min
 Detection: RI
 Temperature: 25 °C
 Samples: Polyethylene Glycol 200 (3 g/L)
 Injection volume: 50 µL

4. Dependence of calibration curves on eluent

4-1. Temperature dependence

Fig. 5 shows the results of an investigation of the temperature-dependence of calibration curves with standard polyethylene oxide, using the TSKgel G5000PW_{XL}-CP, with 0.1 mol/L aqueous solution of sodium nitrate as the eluent. The figure confirms that there is no temperature dependence in the range from 25 °C to 60 °C.

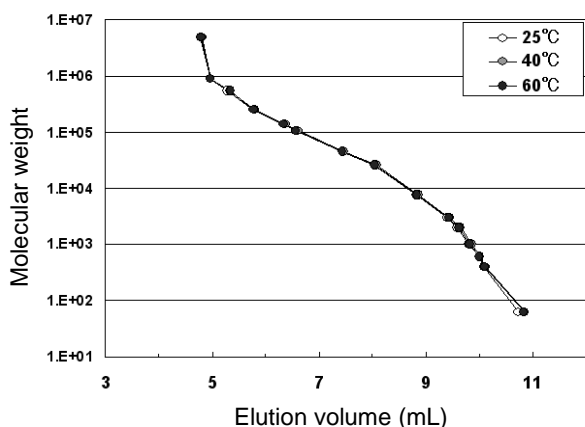


Fig. 5 Dependence of calibration curves on temperature in TSKgel G5000PW_{XL}-CP

Column: TSKgel G5000PW_{XL}-CP (7.8 mm I.D. x 30 cm)
Eluent: 0.1 mol/L sodium nitrate
Flow rate: 1.0 mL/min
Temperature: 25 °C, 40 °C, 60 °C
Samples: Standard polyethylene oxide, polyethylene glycol, ethylene glycol
Injection volume: 100 µL

4-2. Dependence on pH in eluent

Fig. 6 shows calibration curves produced with standard polyethylene oxide when the pH of the eluent (0.1 mol/L acetate buffer) is varied, using the TSKgel G5000PW_{XL}-CP.

The calibration curve showed no changes within a pH range of 4.5 to 8.3.

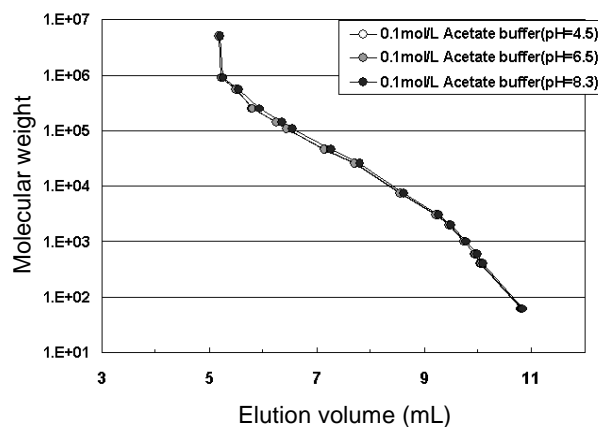


Fig. 6 Dependence of calibration curves on pH of eluent in TSKgel G5000PW_{XL}-CP

Column: TSKgel G5000PW_{XL}-CP (7.8 mm I.D. x 30 cm)
Eluent: 0.1 mol/L acetate buffer (pH = 4.5, 6.5, 8.3)
Flow rate: 1.0 mL/min
Temperature: 25 °C
Samples: Standard polyethylene oxide, polyethylene glycol, ethylene glycol
Injection volume: 100 µL

4-3. Effects of organic solvent concentration in eluent

Fig. 7 shows the results of an investigation of the effect of acetonitrile concentration on the calibration curve using the TSKgel G5000PW_{XL}-CP as the column, when acetonitrile is added to the eluent (0.1 mol/L sodium nitrate). In an acetonitrile concentration range of 0 to 20% very little impact on the calibration curve can be seen. As the acetonitrile concentration increases, elution of the standard samples becomes more rapid, and overall the calibration curve tends to shift toward the low side.

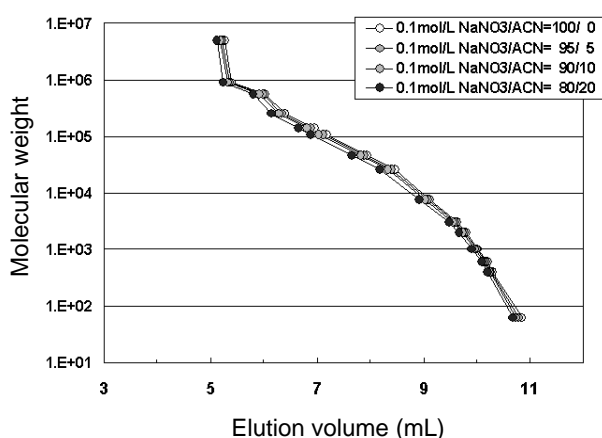


Fig. 7 Effect of concentration of organic solvent on calibration curve of TSKgel G5000PW_{XL}-CP

Column: TSKgel G5000PW_{XL}-CP (7.8 mm I.D. x 30 cm)

Eluent: 0.1 mol/L sodium nitrate/acetonitrile (0/100 to 20/80)

Flow rate: 1.0 mL/min

Temperature: 25 °C

Samples: Standard polyethylene oxide, polyethylene glycol, ethylene glycol

Injection volume: 100 µL

4-4. Flow rate dependency

Fig. 8 shows the results of an investigation of the dependency of standard polyethylene oxide calibration curves on the flow rate, using the TSKgel G6000PW_{XL}-CP column and 0.1 mol/L sodium nitrate as the eluent. The results confirmed that there were no marked changes in the calibration curve over a range of flow rates from 0.25 to 1.0 mL/min.

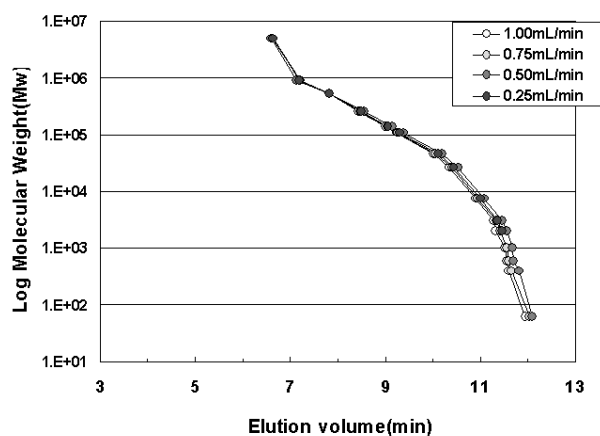


Fig. 8 Dependence of TSKgel G6000PW_{XL}-CP calibration curves on flow rate

Column: TSKgel G6000PW_{XL}-CP (7.8 mm I.D. x 30 cm)

Eluent: 0.1 mol/L sodium nitrate

Flow rate: 0.25 to 1.0 mL/min

Temperature: 25 °C

Samples: Standard polyethylene oxide, polyethylene glycol, ethylene glycol

Injection volume: 100 µL

5. Elution characteristics of cationic polymers

5-1. Comparison of elution using commercial SEC column

Fig. 9 shows the results of a comparison of the elution characteristics of poly(allylamine hydrochloride) (PAA-HCl), a cationic polymer, under mild solvent conditions [neutral solvents with low salt concentration (0.1 mol/L aqueous solution of sodium nitrate)] with the TSKgel G5000PW_{XL}-CP, TSKgel G5000PW_{XL}, and a commercial SEC column.

Due to adsorption, the cationic polymer PAA-HCl completely failed to elute in the PW_{XL} and the commercial aqueous SEC columns. However, the chromatogram obtained with the TSKgel G5000PW_{XL}-CP series showed good elution of the sample.

Fig. 10 shows chromatograms obtained analyzing PAA-HCl with a multi-angle light scattering detector (MALS) using the TSKgel G6000PW_{XL}-CP and TSKgel G6000PW_{XL}, under the same solvent conditions as in **Fig. 9** (0.1 mol/L aqueous solution of sodium nitrate).

The figure verifies that even with MALS, PAA-HCl fails to elute in the TSKgel G6000PW_{XL} due to adsorption and is not detected.

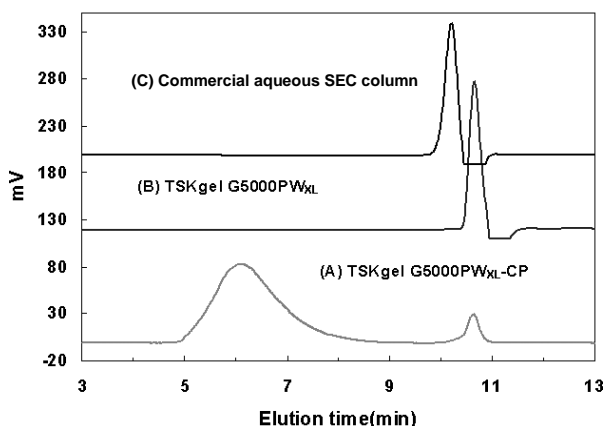


Fig. 9 Chromatograms of PAA-HCl using TSKgel G5000PW_{XL}-CP and conventional columns

Columns: (A) TSKgel G5000PW_{XL}-CP (7.8 mm I.D. x 30 cm)
 (B) TSKgel G5000PW_{XL} (7.8 mm I.D. x 30 cm)
 (C) Commercial aqueous SEC column (7.8 mm I.D. x 30 cm)

Eluent: 0.1 mol/L sodium nitrate

Flow rate: 1.0 mL/min

Detection: RI

Temperature: 25 °C

Samples: Poly(allylamine hydrochloride) (PAA-HCl) (3 g/L)

Injection volume: 100 μ L

Generally, in aqueous SEC, eluents with high salt concentrations are used to inhibit electrostatic interaction with cationic polymers.

Fig. 11 shows a chromatogram produced by repeated analyses of PAA-HCl using a commercial aqueous SEC column in a eluent with a high salt concentration (0.5 mol/L acetate + 0.1 mol/L aqueous solution of sodium nitrate). The figure confirms that although a cationic polymer could be eluted under these solvent conditions, with a small number of injections the sample recovery rate was low, and repeatability was inadequate.

Fig. 12 shows the results of an investigation of the elution characteristics of PAA-HCl using each of the columns of the TSKgel PW_{XL}-CP series, in a 0.1 mol/L aqueous solution of sodium nitrate. Sample recovery rates and average molecular mass (M_w) data obtained by SEC analysis are also shown. The figure confirms that adequate molecular mass data and good recovery rates ($\geq 97\%$) could be obtained with all grades.

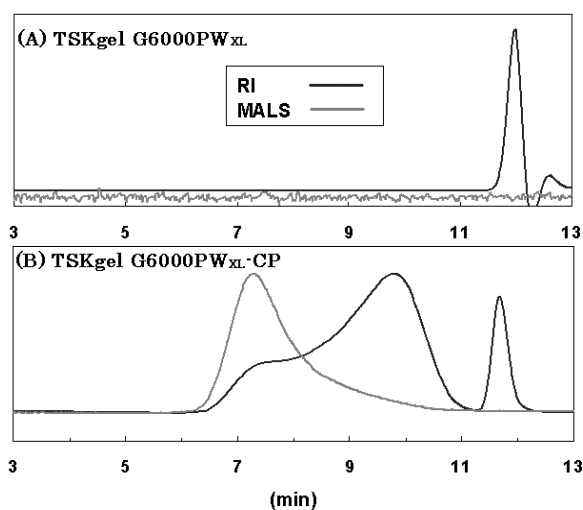


Fig. 10 Chromatograms of PAA-HCl using TSKgel G6000PW_{XL}-CP and TSKgel G6000PW_{XL}

Columns: (A) TSKgel G6000PW_{XL} (7.8 mm I.D. x 30 cm)
 (B) TSKgel G6000PW_{XL}-CP (7.8 mm I.D. x 30 cm)

Eluent: 0.1 mol/L sodium nitrate

Flow rate: 1.0 mL/min

Detection: MALS (DAWN HELEOS), RI

Temperature: 40 °C

Samples: Poly(allylamine hydrochloride) (PAA-HCl) (3 g/L)

Injection volume: 100 μ L

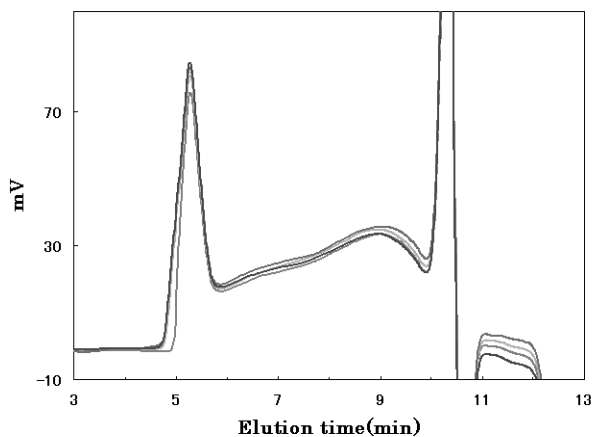


Fig. 11 Chromatograms of PAA-HCl in eluent with high salt concentration using a commercial aqueous SEC column and molecular mass data

Column: Commercial aqueous SEC column
(7.8 mm I.D. x 30 cm)
Eluent: 0.5 mol/L acetic acid + 0.1 mol/L sodium nitrate
Flow rate: 1.0 mL/min
Detection: RI
Temperature: 25 °C
Samples: Poly(allylamine hydrochloride) (PAA-HCl) (3 g/L)
Injection volume: 100 µL

Inject. No.	Mw	Area	Recovery (%)
1	151,000	8,866	91
2	156,000	9,545	98
3	3,300,000	9,650	99
4	3,480,000	9,742	99.6
5	3,510,000	9,778	

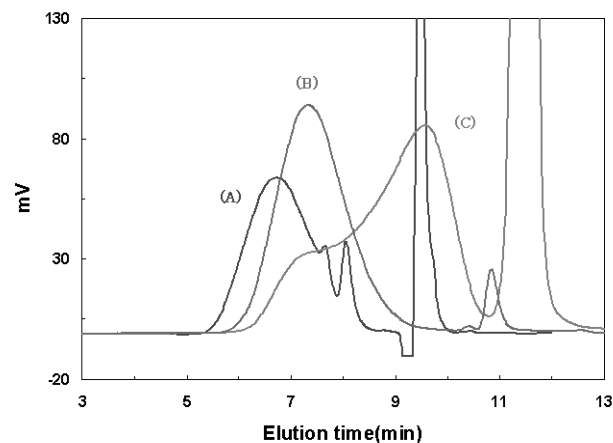


Fig. 12 Chromatograms of PAA-HCl by TSKgel PW_{xL}-CP series

Columns: (A) TSKgel G3000PW_{xL}-CP (7.8 mm I.D. x 30 cm)
(B) TSKgel G5000PW_{xL}-CP (7.8 mm I.D. x 30 cm)
(C) TSKgel G6000PW_{xL}-CP (7.8 mm I.D. x 30 cm)
Eluent: 0.1 mol/L sodium nitrate
Flow rate: 1.0 mL/min
Detection: RI
Temperature: 25 °C
Sample: Poly(allylamine hydrochloride) (PAA-HCl) (3 g/L)
100 µL
(A) PAA-HCl-01
(B) PAA-H-HCl
(C) PAA-HCl-10S

Grade	Mw	Recovery (%)
TSKgel G3000PW _{xL} -CP	6,500	100.2
TSKgel G5000PW _{xL} -CP	168,000	98.8
TSKgel G6000PW _{xL} -CP	455,000	97.4

5-2. Repeatability

Fig. 13 shows the results of an investigation of the repeatability (within a day) of analyses of PAA-HCl in a 0.1 mol/L aqueous solution of sodium nitrate, using the TSKgel G5000PW_{XL}-CP.

From the first run, the sample showed high recovery rates, and the repeatability of chromatograms and average molecular mass (M_w) data analyzed by SEC were confirmed to be excellent.

On the other hand, as shown in **Fig. 14**, repeatability data in the TSKgel G5000PW_{XL} indicated that even with repeated injection of samples under the same conditions, analyte peaks did not appear, clearly indicating that very strong adsorption was occurring.

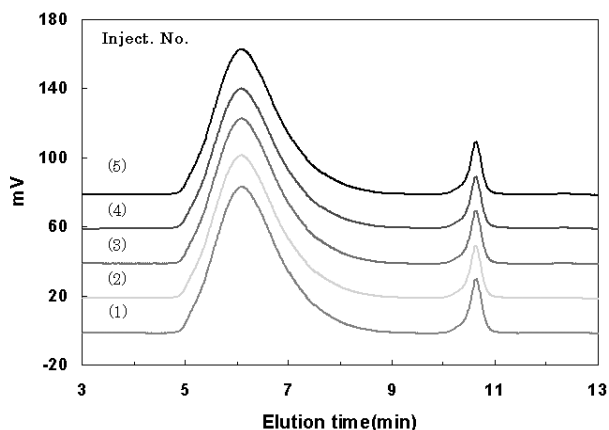


Fig. 13 Repeatability (within a day) of analyses of PAA-HCl using TSKgel G5000PW_{XL}-CP

Column: TSKgel G5000PW_{XL}-CP (7.8 mm I.D. x 30 cm)
 Eluent: 0.1 mol/L sodium nitrate
 Flow rate: 1.0 mL/min
 Detection: RI
 Temperature: 25 °C
 Samples: Poly(allylamine hydrochloride) (PAA-HCl) (0.3%)
 Injection volume: 100 μL

Inject. No.	M _w	Recovery (%)
1	168,000	98.8
2	169,000	99.1
3	168,000	99.1
4	170,000	99.3
5	170,000	99.2

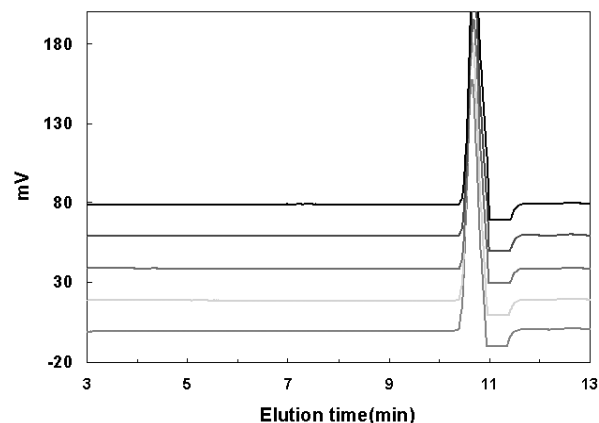


Fig. 14 Repeatability (within a day) of analyses of PAA-HCl using TSKgel G5000PW_{XL}

Columns: TSKgel G5000PW_{XL} (7.8 mm I.D. x 30 cm)
 Eluent: 0.1 mol/L sodium nitrate
 Flow rate: 1.0 mL/min
 Detection: RI
 Temperature: 25 °C
 Sample: Poly(allylamine hydrochloride) (PAA-HCl) (0.3%)
 Injection volume: 100 μL

5-3. Application of various cationic polymers

Fig. 15 shows chromatograms of various polyethyleneimines and poly(allylamine hydrochloride) analyzed with the TSKgel G3000PW_{XL}-CP. **Fig. 16** shows chromatograms of various poly(allylamine hydrochloride) with different molecular masses analyzed with the TSKgel G6000PW_{XL}-CP.

Under each of these conditions, good chromatograms were obtained from the initial run, and samples were eluted in order according to molecular mass (largest molecules first), confirming that SEC elution was occurring normally.

Fig. 17 shows overlapping chromatograms of various cationic polymers with different properties and molecular mass measured in a system in which the TSKgel G6000PW_{XL}-CP, TSKgel G5000PW_{XL}-CP and TSKgel G3000PW_{XL}-CP columns were connected in series.

The TSKgel PW_{XL}-CP series shows an ability to perform good SEC analysis of various cationic polymers with different properties.

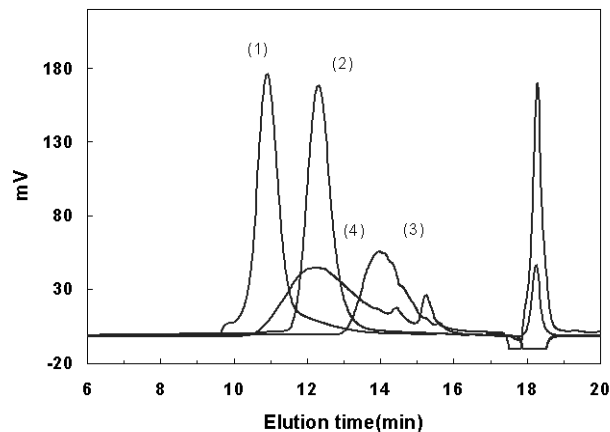


Fig. 15 Chromatograms of various cationic polymers using the TSKgel G3000PW_{XL}-CP

Column: TSKgel G3000PW_{XL}-CP (7.8 mm I.D. x 30 cm x 2)
 Eluent: 0.1 mol/L sodium nitrate
 Flow rate: 1.0 mL/min
 Detection: RI
 Temperature: 25 °C
 Sample: Cationic polymers (3g/L) 100 μ L
 (1) Polyethyleneimine (10000)
 (2) Polyethyleneimine (1800)
 (3) Polyethyleneimine (300)
 (4) Poly(allylamine hydrochloride)

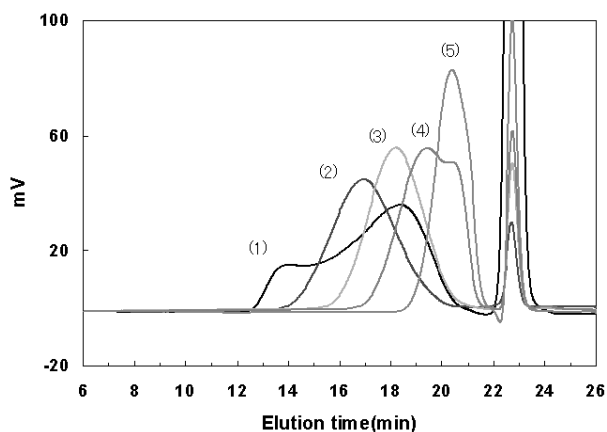


Fig. 16 Chromatograms of various poly(allylamine hydrochloride) using the TSKgel G6000PW_{XL}-CP

Columns: TSKgel G6000PW_{XL}-CP (7.8 mm I.D. x 30 cm x 2)
 Eluent: 0.1 mol/L sodium nitrate
 Flow rate: 1.0 mL/min
 Detection: RI
 Temperature: 25 °C
 Sample: Poly(allylamine hydrochloride) (3g/L) 100 μ L
 (1) PAA-HCl
 (2) PAA-H-HCl
 (3) PAA-HCl-3L
 (4) PAA-HCl-05
 (5) PAA-HCl-01

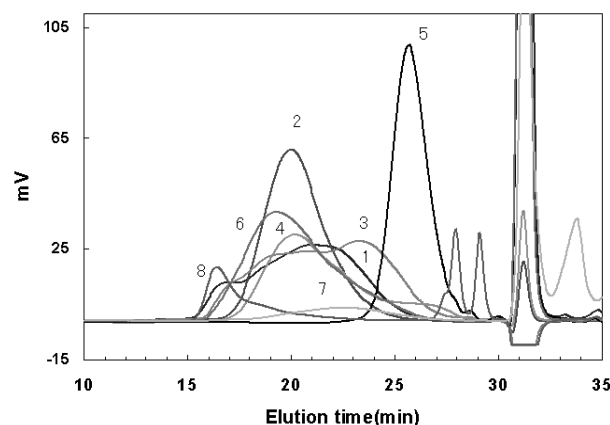


Fig. 17 Chromatograms of various cationic polymers using the TSKgel PW_{XL}-CP series

Columns: TSKgel G(6000 + 5000 + 3000)PW_{XL}-CP (7.8 mm I.D. x 30 cm x 3)
 Eluent: 0.1 mol/L sodium nitrate
 Flow rate: 1.0 mL/min
 Detection: RI
 Temperature: 25 °C
 Sample: Cationic polymers (3g/L) 100 μ L
 (1) Poly(allylamine hydrochloride) (PAA-HCl)
 (2) Poly(allylamine hydrochloride) (PAA-H-10C)
 (3) Polyethyleneimine
 (4) Polydiallyl dimethyl ammonium chloride
 (5) Polydiallyl dimethyl ammonium chloride/sulfur dioxide copolymer
 (6) Polydiallyl dimethyl ammonium chloride/polyacrylamide copolymer
 (7) Cationized dextran
 (8) Chitosan

6. Analysis of absolute molecular mass of cationic polymers by SEC-MALS

The radius of inertia and absolute molecular mass of various cationic polymers were investigated using SEC-MALS (multi-angle light scattering detector) using TSKgel PW_{XL}-CP columns.

6-1. Absolute molecular mass of polydiallyl dimethyl ammonium chloride

Fig. 18 shows absolute molecular mass data (MALS) and a chromatogram (RI) of polydiallyl dimethyl ammonium chloride (PDADM-NH₄Cl) analyzed by SEC-MALS (multi-angle light scattering detector) using the TSKgel G6000PW_{XL}-CP and 0.1 mol/L aqueous solution of sodium nitrate as the eluent.

Fig. 19 shows the relationship between the absolute molecular mass and radius of inertia. It is clear from these results that the polymer components are being eluted in order starting with the largest molecules first, and that good SEC separation is occurring. The relationship between the absolute molecular mass and the radius of inertia is also good.

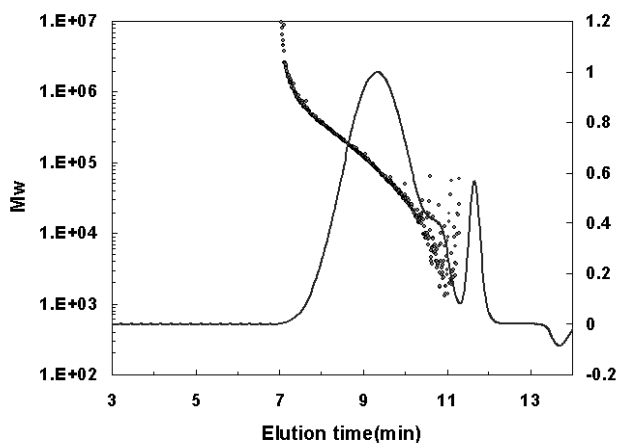


Fig. 18 Absolute molecular mass and chromatogram of elution of PDADM-NH₄Cl by the TSKgel G6000PW_{XL}-CP

Column: TSKgel G6000PW_{XL}-CP (7.8 mm I.D. x 30 cm)
Eluent: 0.1 mol/L sodium nitrate
Flow rate: 1.0 mL/min
Detection: MALS (DAWN HELEOS), RI
Temperature: 40 °C
Sample: Polydiallyl dimethyl ammonium chloride (1 g/L) 100 μ L

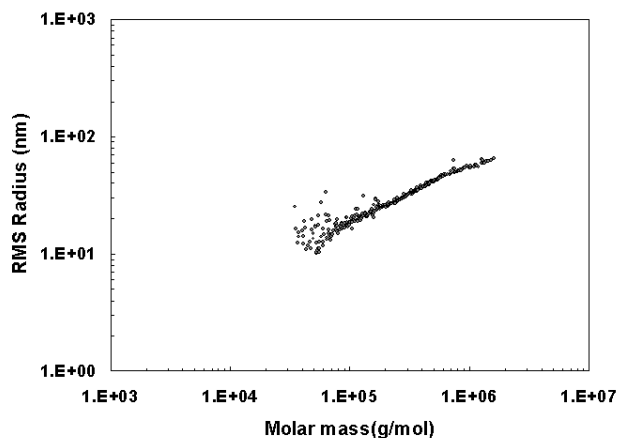


Fig. 19 Relationship between absolute molecular mass and radius of inertia of PDADM-NH₄Cl using the TSKgel G6000PW_{XL}-CP

Same as Fig. 18
Sample: Polydiallyl dimethyl ammonium chloride (1 g/L) 100 μ L

6-2. Absolute molecular mass of PDADM-NH₄Cl copolymer

The absolute molecular mass of copolymers of PDADM-NH₄Cl and acrylamide, and PDADM-NH₄Cl and sulfur dioxide copolymer was analyzed under the same conditions as **Fig. 19**.

Fig. 20 shows the relationship between absolute molecular mass and radius of inertia. It is clear that SEC separation is occurring normally.

As is clear from the data provided here, although the absolute molecular mass of each of these copolymers is approximately the same (about 200,000), the radius of inertia of the copolymer of PDADM-NH₄Cl and acrylamide is about 1.5 times as large.

7. Conclusion

Until now it has been difficult to perform normal SEC analysis of aqueous cationic polymers, because it has been impossible to inhibit electrostatic interaction between the sample and the packing material.

To overcome such problems, we have improved the ionic characteristics on the surface of the packing material using a novel technique, and have confirmed that cationic polymers can now be analyzed with good repeatability, even in neutral eluents with low salt concentrations.

Data from analyses of absolute molecular mass by SEC-MALS also confirms that normal SEC separation is occurring with these columns.

The use of the TSKgel PW_{XL}-CP series is expected to make it easier to analyze various cationic polymers.

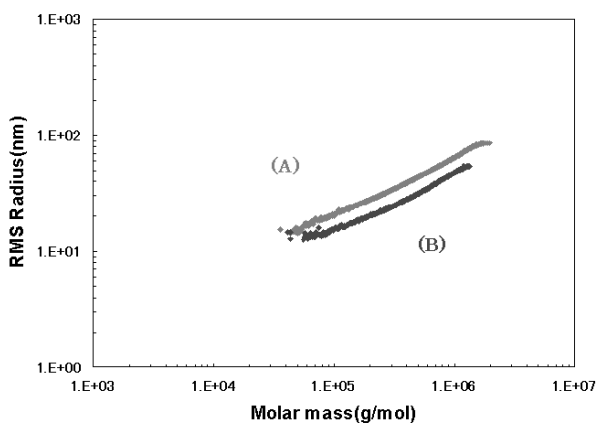


Fig. 20 Relationship between absolute molecular mass and radius of inertia of PDADM-NH₄CL copolymer using the TSKgel G6000PW_{XL}-CP

Same as Fig. 21 [\[Please verify this figure number.\]](#)

Samples: (A) PDADM-NH₄Cl/polyacrylamide copolymer

(1 g/L) 100 μ L

(B) PDADM-NH₄Cl/sulfur dioxide copolymer

(1 g/L) 100 μ L